

Invacare® ProSPIN™ X4™ Ultralightweight Custom Manual Wheelchair



ProSPIN™ X4™

Yes, you can.®

- Forward Design and Breathtaking Responsiveness
- Smart Adjustability



PROSPIN™ X4™

Anything else is a compromise.™

*Invacare® ProSPIN™ X4™
– Folding Performance
Wheelchair.*

*Finally, the next evolution in
performance folding ultra-
lightweight wheelchairs has
arrived. Coupled with an
extremely lightweight frame
and portability, the X4 is the
only choice for someone who
is looking for the chair of the
future. Awe-inspiring design,
premium performance, and
unbelievably lightweight:
The ProSPIN X4 is 360° of
pure exHilaration!*



Future Forward Design and Breathtaking Responsiveness

The ProSPIN X4 takes the look, feel and performance of a folding wheelchair into the new century. From any angle, the X4 appeals to the visual sense like no other wheelchair.

The unique side frame and axle plate design merge aesthetics and performance. The X4 side frame minimizes both welded and bolt-on connections, strengthening the frame and lessening the amount of propulsion energy lost through flex and vibration. The curved axle plate completes the flowing side view, making the X4 look more like a work of art than a wheelchair.

In the front of the frame, the head tube design allows consumers to get closer to desks or other furniture. The caster is positioned closer to the drive

wheel for quickened response during propulsion, and less caster interference with foot placement and maneuverability. The sleek design of the curved wheel lock handle makes it easy to engage and helps maintain a low profile for less interference with propulsion and transfers. Finally, the front capped H-blocks help rigidify the frame by holding the seat rails firmly in place. This greatly reduces the flex associated with many folding-frame wheelchairs.

The ProSPIN X4 is available in a variety of color options. Choose from Contemporary, Too Hot, or Way Cool colors like Tangerine Red and Cosmic Blue.




The unique side frame design puts the X4 in a class by itself.



New swing-in/swing-out front riggings help make transfers easier.



A new curved handle wheel lock adds to the functionality.

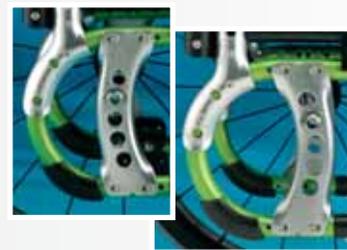

The axleplate can be flipped 180° for plenty of center-of-gravity and seat-to-floor combinations.

Detailed Perfection and Smart Adjustability™

The level of design detail in the available options and accessories further complements the ProSPIN X4 aesthetics. The swing-in/swing-out front riggings allow consumers to swingaway the footrests either to the outside or underneath the wheelchair. This allows easier transfers from the wheelchair into an automobile seat or office chair without having to back up the wheelchair.

The curved axle plate not only looks awesome, it also is incredibly functional, providing plenty of center-of-gravity options and seat-to-floor heights. The axle plate can be moved along the frame in 1/2" increments and has slotted holes to

allow for more center-of-gravity adjustments. If that isn't enough, the entire bracket can be flipped 180° to ensure that the X4 can be configured for optimal stability and performance. An adjustable-angle back option also is available allowing up to 15° of either posterior or anterior adjustment.

With Invacare's own Smart Adjustability™, making adjustments on the fly is surprisingly quick and simple. The adjustable front head tubes allow changes to the front caster angle as needed, without having to remove the front caster wheel. The X4's wheel-locks switch in seconds from push-to-lock to pull-to-lock.

Extremely Lightweight

The ProSPIN X4 was designed to be one of the lightest-weight folding wheelchairs available while also improving the durability needed to endure many years of folding and unfolding. The improved crossbraces allow for a tight, compact fold so the X4 can be easily transported in the back seat or trunk of most vehicles.



Fixed Front Frame Version.

PROSPIN X4

Anything else is a compromise™



New elevating leg rest designed to provide ample adjustability.



The X4 caster head tube adjustment is accessible and simple to change.



A front capped H block helps rigidize the frame by holding the seat rails in place.



The X4 folds tightly for easy storage and transportation.



Yes, you can.®




Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Invacare, the Blue Medallion Design and Yes, you can. are registered trademarks and Anything else is a compromise, ProSPIN, QuickSet, X4 and Smart Adjustability are trademarks of Invacare Corporation. Spinergy is a registered trademark of Spinergy, Inc. Specifications are subject to change without notification.

© 2006
Invacare Corporation
Form No. 05-005
Rev. 09/06

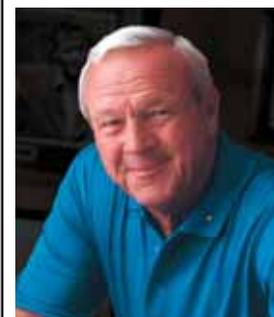

"Invacare says 'Yes, you can.' and I couldn't have said it better myself."

ProSPIN™ Specifications

Frame style	Folding frame with swingaway front
Seat width	14" - 20"
Seat depth	14" - 20"
Back height	10" - 20"
Adjustable-angle back	85°-105° (in 5° increments)
Camber angles	0°, 3°, 6°
Tubing outside diameter	1"
Weight of chair	15 lb.** , 30 lb.***
Weight capacity	250 lb.**
Limited warranty	Limited lifetime on frame and crossbraces

** Base chair weight without accessories.

*** 16" x 16"; AT99, BH20, BRG2, BLUM2, SPCR1, WLPLL, SPIN24, PRIMO, 70STAPER, 1651, HLOOP, IMPACTG, BACK1216, SLIPBACK, 1671LG, ALUMPC, A1577.

Knee-to-heel*

Folding frames w/ seat rails	min	max
	60"	14.00" 19.00"
70"	14.00" 19.00"	
70" Tapered	14.50" 19.50"	

* Range adjusts in 1" increments. Measurement taken from top of seat rail to back of footrest.

Arm heights*

AT99*	8" - 10"
CAD*	10" - 14"
CAF*	10" - 14"
TAD**	6" - 10"
TAF**	8" - 14"

* Range adjusts in 1" increments. ** Range adjusts in 0.5" increments. All measurements from top of seat rail to top of armrest. On AT99, add 1" for upholstered arm pads. For less-seat-rail chairs, add 1.5" to above measurements (top of upper frame rail to top of armrest.)

Rear seat-to-floor heights

13.5" - 18"	20 Com, 20SP
14.5" - 19"	22 Com, 22SP
15.5" - 20"	24 Com7, 24SP, 24HSP, 24SPIN
16.5" - 20.5"	25SPIN, 559TS, 26SP, 26HSP, 26SPIN

Front Seat-to-floor heights*

Fork	Spacer	3"	4"	5"	6"
			4SPLT AT900	5SOFT AT901 4ALUM2	AT902 5ALUM2
BRG1	SPCR0	15.50" 16.00"	15.50" 16.00" 16.50"	16.50" 17.00"	17.50"
	SPCR1	16.50" 17.00"	16.50" 17.50"	17.50" 18.00"	18.50"
	SPCR2	17.50" 18.00"	17.50" 18.50"	18.50" 19.00"	19.50"
BRG2	SPCR0	17.50" 18.00"	17.50" 18.00" 18.50"	18.00" 18.50" 19.00"	18.50" 19.00" 19.50"
	SPCR1	18.50" 19.00"	18.50" 19.50"	19.00" 20.00"	19.50" 20.50"
	SPCR2	19.50" 20.00"	19.50" 20.50"	20.00" 21.00"	20.50" 21.50"
SOFR-F	SPCR0	16.50" 17.00"	16.50" 17.50"	17.50" 18.00"	18.50"
	SPCR1	17.50" 18.00"	17.50" 18.50"	18.50" 19.00"	19.50"
	SPCR2	n/a	n/a	n/a	n/a