




MaTRx seating series

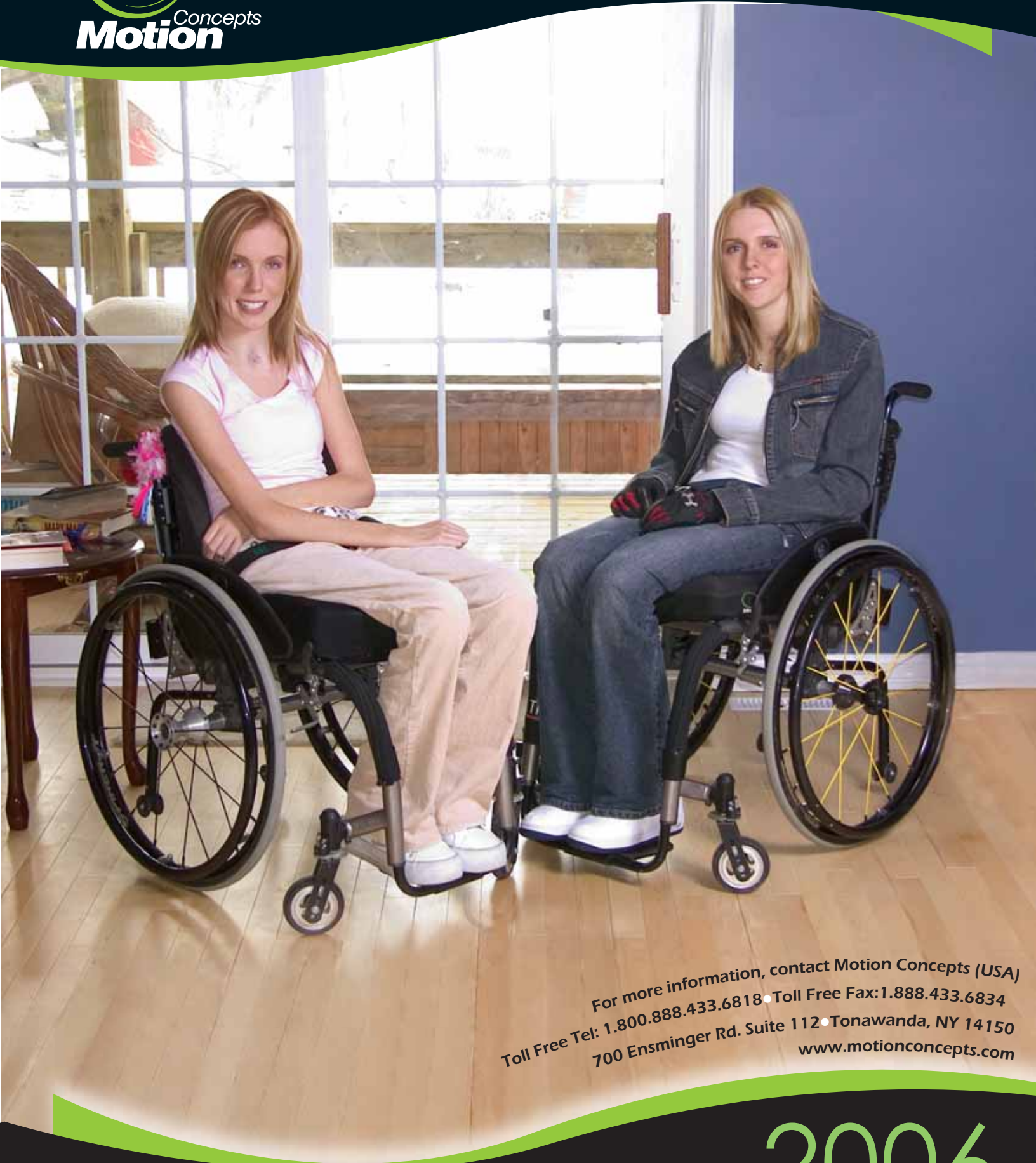

For more information, contact Motion Concepts (USA)
Toll Free Tel: 1.800.888.433.6818 • Toll Free Fax: 1.888.433.6834
700 Ensminger Rd. Suite 112 • Tonawanda, NY 14150
www.motionconcepts.com

2006

NEW!

MULTI-POSITION ARM PADS



- ▶ Soft, gel-like feel provides extra comfort and protection for arms
- ▶ Highly durable skinned polyurethane resists cracking and peeling
- ▶ Triple track extrusion allows adjustability of mounting position and attachment location for other accessories
- ▶ 3" wide; desk length (10") or full length (14")

ERGONOMIC ARM TROUGH/ HAND PAD KIT



- ▶ Molded, skinned polyurethane fits the natural arm curves length
- ▶ Unique cut-out at elbow provides users with flexibility of various arm positions
- ▶ Triple track extrusion allows a variety of mounting positions
- ▶ Complete with flat hand pad or interfaces with other manufacturer's hand pieces

FOAM-IN-PLACE KIT



- ▶ Easy to use chemical kit allows quick customization of MaTRx PB Backs to accommodate various user needs
- ▶ Individually molded foam insert produced in minutes with user seated in back

NEW!

MATRIX-V HEAVY DUTY CUSHION



- ▶ Weight Capacity 650 lbs
- ▶ Widths up to 30", depths to 26"
- ▶ Cushion height (at leg trough) 2 1/4"



PELVIC OBLIQUITY KIT (POK)

- ▶ Easily installed in all MaTRx-V and Kid*ab*ra cushions/bases
- ▶ Accommodates 1/2" or 1" obliquities



POSTURE BACK CANE REPLACEMENT

- ▶ Replaces existing back canes which are less than 1" diameter for use with any MaTRx Backs




EXTENDED SHEAR REDUCTION (ESR) MOUNTING HARDWARE

- ▶ Allows the MaTRx PB Back to be mounted directly to Motion Concepts TRx power seating with ESR



SEAT DEPTH REDUCTION/ LEG LENGTH DISCREPANCY

- ▶ Foam, inner and outer covers modified to your specifications upon request

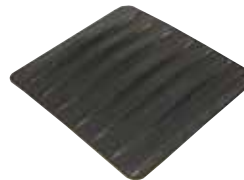

SEAT RIGIDIZER

NEW!

- ▶ Molded plastic solid seat insert
- ▶ Eliminates hammock effect of sling upholstery to provide stable support
- ▶ Lightweight, slim profile
- ▶ Adheres securely to cushion inside cover
- ▶ Sizes 14" - 20" W x 16" - 20" D



DURA CLIP

NEW!

- ▶ Mounting clips provide secure lock over drop hooks
- ▶ Easy lock and release button
- ▶ Highly durable
- ▶ Constructed from bullet-proof composite material

matrx flo-tech CUSHION

Cushions Shaped Around People

NEW!

- ▶ Market leading fluid pressure reduction cushion designed for users with highest risk of skin breakdown or existing tissue damage
- ▶ Unique dual layer fluid sacs provide immersion to sensitive bony prominences and prevent "bottoming out"
- ▶ Comfortable, molded foam base includes pre-ischial contour and trochanteric shelf
- ▶ Moisture resistant "Platilon" inner liner protects foam from moisture
- ▶ Stretchy, water resistant, breathable outer cover – includes hook and loop fasteners, lifting strap and side pocket
- ▶ Optional fluid pack supplements/add-on components available to accommodate individual postural needs
- ▶ Sizes from 15" - 20" W x 16" - 21" D; 300 lb weight capacity



matrx-v CUSHION/BASE

Effective, Reliable Pressure Reduction



- ▶ Excellent pressure reduction, comfort and positioning designed for users at moderate to high risk of skin breakdown
- ▶ Anatomically designed molded Visco-elastic/HR foam combination includes pre-ischial contour, trochanteric shelf, beveled front-end
- ▶ Moisture resistant, zippered inner cover protects foam and allows easy access for cushion modification
- ▶ Moisture resistant/breathable outer cover available in black or appealing Kid*ab*ra fabric – includes hook and loop fasteners, lifting strap, side pocket
- ▶ Optional Pelvic Obliquity Kit/ Leg Length Discrepancy/Seat Depth Reduction
- ▶ Lightweight, maintenance free
- ▶ Standard sizes from 12" - 22" W x 12" - 20" D, 300 lb weight capacity
- ▶ New Heavy Duty sizes available – up to 30" W and 26" D; 650 lb weight capacity

MaTRx-V
Cushion/Base



Kid*ab*ra
Cushion/Base



NEW!

matrx-pb elite

POSTURE BACK

Support for Active Lifestyles

- ▶ "Easy Set" hardware allows ease of adjustment and optimal orientation of back in wheelchair
- ▶ Very lightweight – most models under 3 lbs!
- ▶ Sleek, stylish and breathable back shell with cut-outs for breathability
- ▶ 3" contour depth suits more active users
- ▶ Wide range of sizes from 12" to 20" W x 12", 16" or 20" H
- ▶ Optional formable, width/angle adjustable lateral supports
- ▶ 300 lb weight capacity



'Easy Set' hardware used for both the MaTRx-PB Elite and Posture Back

matrx-pb

POSTURE BACK

Infinitely Adjustable Postural Control

- ▶ "Easy Set" hardware installed and adjusted in minutes – while user is seated in wheelchair
- ▶ Infinite adjustment of recline (30°), seat depth (4") and even rotation!
- ▶ Each back fits 3 different wheelchair widths
- ▶ 16" and 20" high backs include mounting plate/attachment point for headrests
- ▶ Heights 12", 16" or 20"; widths from 14" - 24"
- ▶ Available in Standard 5" or Deep 7" contour depth
- ▶ Optional Visco Foam upgrade or Foam-in-Place kit
- ▶ Weight capacity: 300 lbs; Heavy Duty: 500 lbs



matrx

CONTOUR BACK

NEW!

- ▶ Simple and comfortable basic back support
- ▶ Mildly contoured, limited depth/recline adjustability
- ▶ Pre-drilled for attachment of headrest supports
- ▶ Fixed hardware mounting ensures back remains securely attached to wheelchair
- ▶ Sizes: 15" - 20" W x 12", 15" or 17" H
- ▶ Weight capacity: 300 lbs

