


# Jazzy® 1170XL

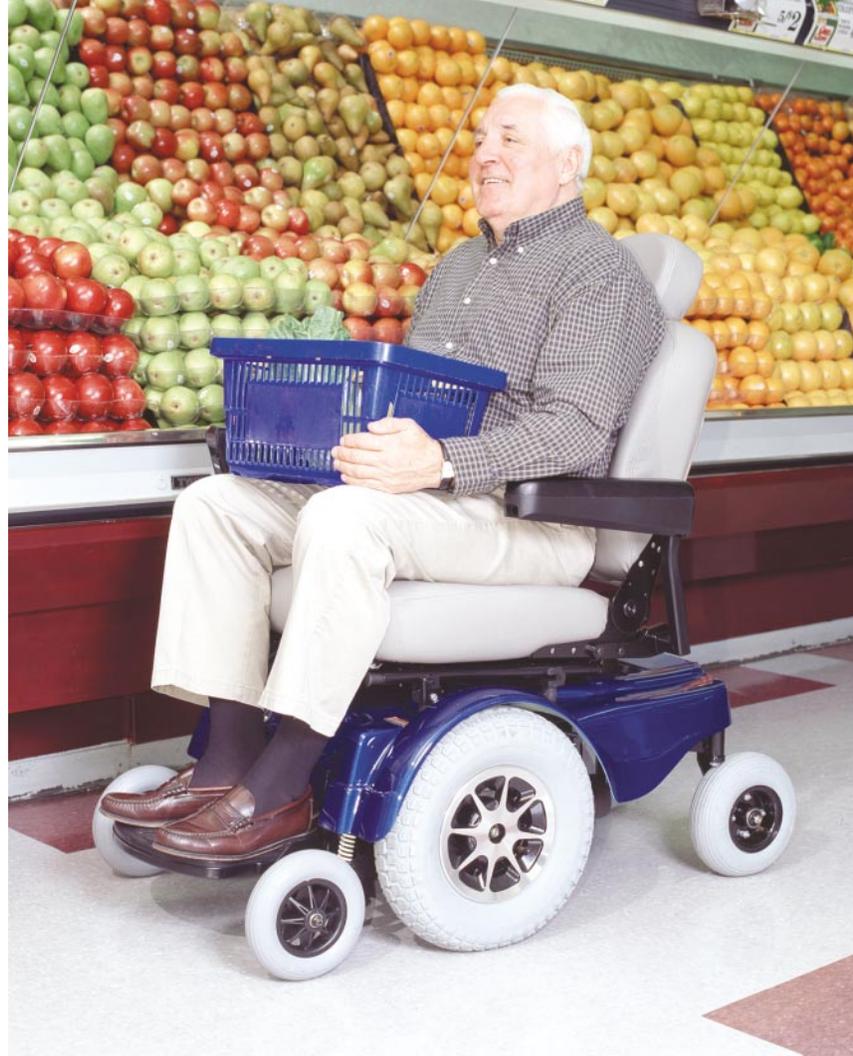

Jazzy 1170 XL shown in Candy Apple Red with optional high-back seat, standard PG Pilot controller and foot platform

The Jazzy® 1170 XL is a stylish, high performance power chair that combines excellent indoor maneuverability with larger motors and a longer frame to allow for a 30% increase in torque, a top speed of 6 mph and a 400 lbs. weight capacity. The Jazzy 1170 XL comes standard with


Active-Trac® Suspension which provides for great outdoor performance and rough terrain handling. The Jazzy 1170 XL is the ultimate power chair for the active outdoor user who still requires tight-quarter maneuverability.

**Pride**®  
Mobility Products Corp.

(US) 800-800-8586 · (Canada) 888-570-1113 · [pridemobility.com](http://pridemobility.com)

<b>Weight capacity:</b>	400 lbs.
<b>Maximum speed:</b>	Up to 6 mph
<b>Ground clearance:</b>	4"
<b>Turning radius:</b>	23"
<b>Overall length:</b>	46"
<b>Overall width:</b>	26.5"
<b>Wheels:</b>	Drive: 16" pneumatic Casters: 9" pneumatic or solid Anti-tips: 8" solid
<b>Suspension:</b>	Active-Trac®
<b>Drive train:</b>	Two-motor, mid-wheel drive
<b>Braking system:</b>	Intelligent braking (electronic, regenerative disc brakes)
<b>Std. Electronics:</b>	70 amp, PG Pilot controller
<b>Specialty controls:</b>	Yes
<b>Battery charger:</b>	On-board, 5 amp (standard)
<b>Per-charge range:</b>	Up to 25 miles
<b>Battery requirements:</b>	(2) 12 volt, deep cycle
<b>Battery size:</b>	70AH, Group 24
<b>Transport weight</b>	
<b>Battery weight:</b>	53.5 lbs. each
<b>Base weight:</b>	132 lbs.
<b>Std. seat weight:</b>	33 lbs. (medium-back)
<b>Seat dimensions:</b>	24" x 26" (max. contoured seating)
<b>Warranty:</b>	Lifetime limited warranty on frame 2-year limited warranty on electronics 18-months limited warranty on drive motors

**Colors:**

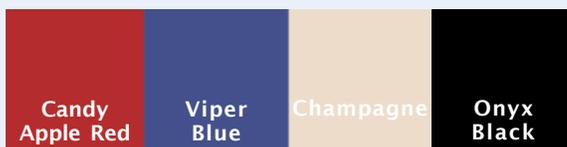

The information contained herein is correct at the time of publication; We reserve the right to alter specifications without prior notice. Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure.

FEATURES

- Powerful, patented, two-motor mid-wheel drive design provides superb power and acceleration
- Active-Trac® and rear caster suspension for outstanding outdoor performance
- Adjustable front anti-tips to meet your preferences and terrain challenges
- On-board battery charger
- Battery-conserving sleep mode available with Remote Plus controller

OPTIONS

- Synergy™ seating systems
- TRU-Balance™ Power Positioning Systems
- Solid seat pan
- Power elevating seat
- Flat-free tires
- Lap belt
- Cup holder
- Heavy-duty dust cover
- Swing-away joystick
- Swing-away footrests
- Elevating leg rests
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Rear basket
- Specialty controls
- Angle-adjustable footplates
- 70° and 80° heavy-duty drop-in footrests
- Saddlebag

**1170 XL PLUS PACKAGE:**

- 16" flat-free drive wheels
- 9" solid caster wheels
- 70AH, Pilot controller (standard)
- 80AH, PG Remote Plus or Dynamic controller (optional)
- Contoured high-back or medium-back seat available up to 28" (width) by 26" (depth)
- Variable speed: up to 4.5 mph
- Overall width: 30.5"
- Weight capacity up to 450 lbs.
- Base weight of 136.5 lbs.