

JET™ 2



Jet™ 2

Shown in Jet Red with gray medium-back seat with foot platform and Dynamic controller.

Style, performance and value unite to make the Jet™ 2 one truly impressive power chair. Pride's patented JZMD® mid-wheel drive technology and sharp style bring the Jet 2 the performance and eye-catching look you demand. And the Jet 2's attractive price brings you the value you deserve.

Pride
Mobility Products Corp.

800-800-8586 (Exeter, PA) • 888-570-1113 (Canada) • www.pridemobility.com

FEATURES


JET 2

- Available in Jet Red and Jet Blue
- Powerful two-motor, mid-wheel drive design
- Patented mid-wheel drive design offers outstanding maneuverability
- 14" knobby tires and optimum weight distribution provide superior traction
- Low center of gravity makes handling smooth and stable
- Speed range adjustable up to 4.5 mph
- Adjustable seat height and back angle
- On-board battery charger

JET 2 WITH HD PACKAGE

- PG VSI 50 amp controller (with standard 70 amp boost that engages automatically when needed)
- Available contoured seating up to 22" x 20"
- Up to 350 lbs. weight capacity
- 14" solid drive wheels (required)
- Base weight 97.5 lbs.
- Width 24.5"

COLOR OPTIONS

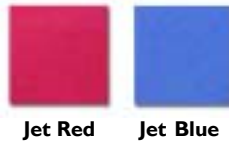

SPECIFICATIONS

Drive Wheels:	14" pneumatic or solid
Rear Casters:	8" solid
Anti-tips:	6" solid, front-mounted
Suspension:	Limited
Maximum Speed:	Up to 4.5 mph
Brakes:	Intelligent Braking (electronic, regenerative disc brakes)
Ground Clearance:	3"
Turning Radius:	19.5"
Overall Size:	Length: 38" Width: 23.5"
Seating Options:	Contoured medium-back, gray or beige vinyl (standard)
Drive Train:	Two-motor, mid-wheel drive
Battery:	Two 12 volt, NF-22
Per-charge Range:	Up to 25 miles
Battery Charger:	4 amp, on-board
Electronics:	Dynamic DL 50 amp
Weight Capacity:	300 lbs.
Weight:	Base: 96 lbs. Seat: 36.5 lbs.
Colors:	Jet Red and Jet Blue
Warranty:	5-year limited warranty on frame 1-year limited warranty on electronics 1-year limited warranty on drive motors



OPTIONS

- Contoured high-back seat, gray or beige vinyl
- Lap belt
- Cup holder
- Heavy-duty dust cover
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Swing-away footrests
- Elevating leg rests
- Heel loops
- Height-adjustable armrests

800-800-8586 (Exeter, PA) • 888-570-1113 (Canada) • www.pridemobility.com

Pride
Mobility Products Corp.