

# QUANTUM<sup>®</sup>

## 1107



The Quantum<sup>®</sup> 1107 easily disassembles to travel with you wherever you go. Yet, complete with advanced Active-Trac<sup>®</sup> Suspension and powerful new GLIDE motors for added power and torque, the Quantum 1107 offers the stability and performance of a full-size, rigid frame chair.

QUANTUM<sup>®</sup> 1107 with  
TRU-Balance<sup>®</sup> Power  
Tilt System



# Q1107

SPECIFICATIONS	
<b>DRIVE WHEELS:</b>	10" pneumatic
<b>REAR CASTERS:</b>	6" solid
<b>ANTI-TIPS:</b>	5" solid, adjustable, front-mounted
<b>SUSPENSION:</b>	Advanced Active-Trac®
<b>MAXIMUM SPEED:</b>	Up to 4.5 mph <sup>1</sup>
<b>BRAKES:</b>	Intelligent Braking (electronic, regenerative disc brakes)
<b>GROUND CLEARANCE:</b>	2.5"
<b>TURNING RADIUS:</b>	20.5"
<b>OVERALL SIZE:</b>	Length: 36.5" Width: 23"
<b>SEATING:</b>	Synergy® Seat Adjustable Width: 10"-20" Adjustable Depth: 10"-20" TRU-Balance® Power Positioning System <sup>2</sup> Adjustable Width: 14"-20" Adjustable Depth: 14"-20" Specialty Seat
<b>SPECIALTY CONTROLS:</b>	Yes
<b>DRIVE TRAIN:</b>	Two-motor, mid-wheel drive
<b>BATTERIES:</b>	Two 12 volt deep cycle U-1
<b>RANGE PER CHARGE:</b>	Up to 25 miles <sup>1</sup>
<b>BATTERY CHARGER:</b>	Off-board, 5 amp (standard)
<b>ELECTRONICS:</b>	75A, Pride Flight
<b>WEIGHT CAPACITY:</b>	300 lbs. maximum (user weight)
<b>BASE WEIGHT:</b>	112 lbs.
<b>SEAT WEIGHT:</b>	37 lbs. Specialty Seat
<b>BATTERY WEIGHT:</b>	24.5 lbs. (each)
<b>COLOR CHOICES:</b>	Candy Apple Red, Viper Blue, Onyx Black, and Sunburst Orange
<b>WARRANTY:</b>	Lifetime limited warranty on frame 2-year warranty on electronics 18-month warranty on drive motor

1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.  
2) 200 lbs. weight capacity



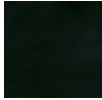

FEATURES
<ul style="list-style-type: none"> <li>• GLIDE motors for added power and torque</li> <li>• Disassembles for ease of transportation</li> <li>• Aluminum frame rails and center locking mechanism provide a very rigid platform for ultimate stability</li> <li>• Innovative battery boxes and front module eliminate battery and motor harnesses</li> <li>• Advanced Active-Trac® Suspension for enhanced performance on rough and uneven terrain, creating a smooth and stable ride</li> <li>• 75 amp, Pride Flight controller</li> <li>• Standard 5 amp, Power Cube off-board charger</li> </ul>

OPTIONS	
<ul style="list-style-type: none"> <li>• Contoured high-back seat</li> <li>• Solid tires</li> <li>• Lap belt</li> <li>• Cup holder</li> <li>• Heavy-duty dust cover</li> <li>• Swing-away joystick</li> <li>• Elevating leg rests</li> </ul>	<ul style="list-style-type: none"> <li>• Oxygen tank holder</li> <li>• Cane/crutch holder</li> <li>• Walker holder</li> <li>• Rear basket</li> <li>• Angle-adjustable footplates</li> <li>• 70° and 80° heavy-duty drop-in footrests</li> </ul>

## QUANTUM 1107

Shown in Candy Apple Red with TRU-Balance® Power Tilt Seating System and Remote Plus controller



AVAILABLE COLORS			
			
Candy Apple Red	Viper Blue	Onyx Black	Sunburst Orange



888-570-1113 • St. Catharines, Ontario • [quantumrehab.com](http://quantumrehab.com)

It is recommended that you consult a qualified Rehab Provider prior to fitting or ordering any Quantum Rehab product. The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.