

QUANTUM[®]

600



EQUIPPED WITH
MID-WHEEL 6[™]
TECHNOLOGY

Equipped with our patent-pending Mid-Wheel 6 Technology, the Quantum 600 delivers maximum stability without compromising performance. The independent active OMNI-Casters (nylon spherical-shaped casters) raise and lower to compensate instantly for changes in the terrain and allow for smooth and stable transitions. In addition, Quantum Rehab has equipped the unit with powerful new GLIDE Motors for added power and torque.




600

SPECIFICATIONS

DRIVE WHEELS:	14" pneumatic
REAR CASTERS:	5" solid with wheel sphere
ANTI-TIPS:	6" solid with wheel sphere
SUSPENSION:	ATX Suspension (Active-Trac® with extra stability)
MAXIMUM SPEED:	Up to 5 mph ¹
BRAKES:	Intelligent Braking (electronic, regenerative disc brakes)
GROUND CLEARANCE:	3"
TURNING RADIUS:	19.5"
OVERALL SIZE:	Length: 35.25" Width: 23.5"
SEATING:	Synergy® Seat Adjustable Width: 10"-20" Adjustable Depth: 10"-20" TRU-Balance® Power Positioning System Adjustable Width: 14"-20" Adjustable Depth: 14"-20" Specialty Seat Kit
DRIVE TRAIN:	Two-motor, mid-wheel drive
BATTERIES:	Two 12 volt deep cycle NF-22
RANGE PER CHARGE:	Up to 25 miles ¹
BATTERY CHARGER:	Lester 8A, Off-board (standard)
ELECTRONICS:	75A, Flight controller
WEIGHT CAPACITY:	300 lbs. maximum (user weight)
BASE WEIGHT:	124 lbs.
BATTERY WEIGHT:	38 lbs. (each)
COLOR CHOICES:	Candy Apple Red, Viper Blue, Sunburst Orange Onyx Black and Champagne
WARRANTY:	Lifetime limited warranty on frame 2-year warranty on electronics 18-month warranty on drive motor

¹) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.

AVAILABLE COLORS

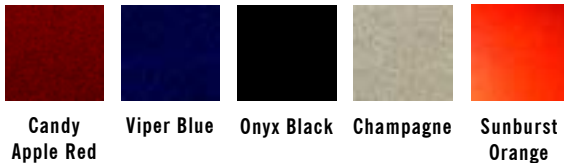

FEATURES

- Mid-wheel 6 Technology allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® Power Positioning System
- OMNI-Casters (nylon spherical-shaped casters) on front and rear prevent wheel hang-ups
- Side-mounted, easily accessible freewheel levers
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Casters for enhanced performance over more varied terrain
- Easy front access to batteries
- Easily accessible side charger port
- New 75 amp, Flight controller or 70 amp PG VSI Integral controller
- Powerful new GLIDE Motors for added power and torque
- Lester 8 amp, off-board charger

OPTIONS

- Solid tires
- Swing-away joystick
- Multi-axis joystick mount
- Elevating leg rests
- Oxygen tank holder
- Neoprene leg rest sleeves
- Gel abductor buttons
- Angle-adjustable footplates
- 70° drop-in footrests
- Lap belt
- On/off-board charger combo

QUANTUM 600

Shown in Candy Apple Red with TRU-Balance® Power Tilt Seating System and Remote Plus controller

