



QUANTUMTM DYNAMO[®]

THE POWER CHAIR FOR PLAYING



QUANTUMTM DYNAMO[®]

THE POWER CHAIR FOR PLAYING

	DYNAMO [®]	DYNAMO [®] ATS
Wheels	10" pneumatic	10" pneumatic
Rear Casters	8" solid	8" solid
Anti-tips	6" solid, front-mounted	6" solid, front-mounted
Suspension	Limited	Active-Trac
Maximum Speed	Up to 4 mph	Up to 4 mph
Brakes	Intelligent Braking	Intelligent Braking
Ground Clearance	2"	2.5"
Turning Radius	16.5"	18"
Overall Size	Length: 32" Width: 23"	Length: 36.5" Width: 23"
Drive Train	Two-motor, mid-wheel drive	Two-motor, mid-wheel drive
Battery	Two 12-volt, U-1	Two 12-volt, U-1
Per-charge Range	Up to 25 miles	Up to 25 miles
Battery Charger	On-board	On-board
Electronics	PG VSI	PG Remote Plus
Weight Capacity	160 lbs.	160 lbs.
Transport Weight	Base: 79.5 lbs. Versa Seat: 44 lbs.	Base: 90.5 lbs. Versa Seat: 44 lbs.

It is recommended that you consult a qualified Rehab Provider prior to fitting or ordering any Quantum Rehab product.



www.quantumrehab.com

Pride

Mobility Products Co.

380 Vansickle Road
St. Catharines, Ontario L2R 6P7
1-888-570-1113

QUANTUM[®]
REHAB
Moving at the Speed of Life[®]
A Division of Pride Mobility Products[®] Corp.

A Division of Pride Mobility Products[®] Corp.



QUANTUM DYNAMO

The Power Chair for Playing

Quantum™ Rehab realizes the potential of the versatile Quantum product line to serve the needs of children with a variety of mobility issues. The new Quantum Dynamo® and Dynamo ATS are the first pediatric power bases to offer Pride's patented mid-wheel drive

design. Both the Dynamo and Dynamo ATS are available with the pediatric Synergy™ seating system, fully programmable electronics, pediatric swing-away footrests and the dependability of a Quantum product.

Why choose the Quantum™ Dynamo® for pediatrics?

The Quantum Dynamo is an ideal power base for pediatric use. Its compact design fits easily into the close quarters of a child's environment, while its 16.5" turning radius gives the unit smooth maneuverability

in the home, at school or on the playground. To keep children eye-to-eye with their friends, both the Dynamo and Dynamo ATS offer an optional seat-to-floor height as low as 15.5".

The Quantum™ Dynamo® ATS

The Quantum Dynamo ATS offers all of the same features as the Dynamo, with the added benefit of Active-Trac Suspension, while maintaining an impressive 18" turning radius. Pride's Active-Trac Suspension allows active children the opportunity to negotiate the roughest terrain to effortlessly keep up with their friends and family.

Quantum Dynamo ATS with Snug Seat Panda in Tilt/Recline Position




Quantum Rehab is devoted solely to customizing the Quantum Power Chair line with specialty seating systems, controllers and a variety of rehab accessories based on an individual's needs and preferences. Our goal is to see that no child's needs go unmet.

Quantum Dynamo ATS shown in Gator Green with Synergy™ seating system and PG Remote Plus controller



Other features of the Dynamo and Dynamo ATS include:

- Stylish and sporty non-medical appearance.
- Seating systems as small as 10" x 10" that can grow with the child.
- Available with all of Quantum Rehab's specialty drive controls.
- Able to interface with a variety of power positioning systems including tilt, tilt/recline and tilt/recline with elevating seat combos.
- 10" knobby drive wheels and patented mid-wheel drive technology offer superior traction over most surfaces that today's active child wants to explore.
- The Dynamo and Dynamo ATS come in a variety of fun colors.

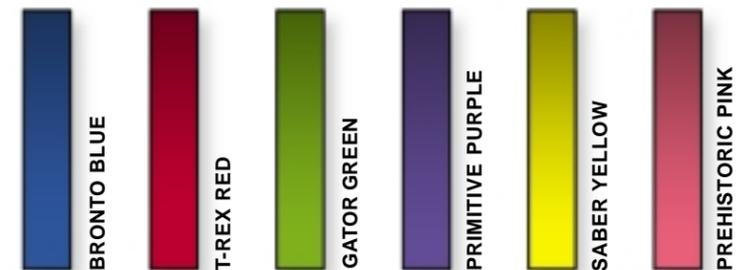

A special Colorbook White shroud that can be colored with dry-erase markers is also available.

*Colors may vary slightly from those shown.



Tiger Folding Mobility Frame
The frame folds for ease of storage and transport.