

INVACARE[®] CANADA SCOOTERS



Yes, you can.™



Yes, you can.®

Invacare® scooters are designed by a team of engineers and clinicians, with input from casual and everyday riders. The result, is a family of scooters that check off every box in terms of performance, design and comfort.

Featuring everything from indoor micro-portable scooters, to outdoor robust, heavy duty scooters. The Invacare Canada scooter lineup is unparalleled in its offering and presents a solution for every need.

Tested at the world class TÜV certification facility in Germany, Invacare scooters have been certified for safety, function and quality for peace-of-mind reliability and long-lasting durability.



Invacare® Canada Scooters

Colibri™ 3-Wheel & 4-Wheel Micro-Portable

A compact scooter ideal for indoor spaces, trips to the mall, visits to the park, or weekends away.

Leo™ 3-Wheel & 4-Wheel Mid-Sized

The perfect indoor/outdoor scooter. Compact and nimble enough to navigate indoor spaces, while offering the range and comfort benefits of a larger scooter.

Pegasus^{METRO} 3-Wheel & 4-Wheel

Ergonomically-superior outdoor model equipped with advanced safety features.

Pegasus^{PRO} 4-Wheel

Elevated design features including an information-rich LCD screen, enhanced suspension and room to accommodate larger batteries for extended travel.

Comet^{PRO} 4-Wheel

A unique suspension system and larger tires make it ideal for tackling rougher outdoor travel terrain.

Comet^{ULTRA} 4-Wheel

Robust model with an increased weight capacity and 26" wide seat with pneumatic lumbar adjustment.





reddot award 2014
winner

COLIBRI4B-12

Colibri™ 3-Wheel & 4-Wheel Micro-Portable

Model Numbers:

	3-Wheel 12 Ah Batteries	4-Wheel 12 Ah Batteries	4-Wheel 18 Ah Batteries
Orange	COLIBRI3O-12	COLIBRI4O-12	COLIBRI4O-18
Red	COLIBRI3R-12	COLIBRI4R-12	COLIBRI4R-18
White	COLIBRI3W-12	COLIBRI4W-12	COLIBRI4W-18
Green	COLIBRI3G-12	COLIBRI4G-12	COLIBRI4G-18
Blue	COLIBRI3B-12	COLIBRI4B-12	COLIBRI4B-18
Silver	COLIBRI3S-12	COLIBRI4S-12	COLIBRI4S-18

The Invacare® Colibri is the perfect travel companion. Designed with unparalleled durability features and enhanced with patented LiteLock™ Technology for simple disassembly, the Colibri is the ideal solution for trips to the mall, visits to the park or weekends away. Standard with batteries and a 300 lbs. weight capacity.

Features

- Three platforms with different wheel and battery options: 200 x 50 three-wheel with 12 Ah batteries, 200 x 50 four-wheel with 12 Ah batteries or 210 x 65 four-wheel with 18 Ah batteries.
- Patented LiteLock™ technology enables easy assembly and disassembly for transportation and travel.
- Six interchangeable colour panel options to personalize the Colibri to suit your unique style.
- Durable, double-spoke black wheels and composite front basket.
- Adjustable three-position tiller with delta steering handle.
- Padded vinyl driver seat featuring width and angle adjustable arms, height adjustment and stylish accent stitching.
- 12 Ah and 18 Ah dual battery pack options offer convenient charging through the front port and up to a 16 km travel range.
- Provider-programmable electronics to reverse throttle direction and adjust driving performance.
- Standard built-in safety features include rear anti-tippers, front and rear mounted reflectors and brake disengagement lever.
- 300 lbs. weight capacity.

STYLE IS FOREVER. INVACARE SCOOTERS.

Colibri™ Scooter Features



Six Colour Panel Options

Personalize your scooter with one of six interchangeable colour options available. Choose from Diamond Silver, Sunstone Orange, Jade Green, Jasper Red, Sapphire Blue or Marble White.



LiteLock™ Technology

Colibri's patented LiteLock system allows you to simply remove the battery box to disengage the locking mechanism for easy disassembly of the scooter into manageable components for travel.



Comfortable and Adaptable

The three-position tiller with delta steering handles ensure a comfortable driving position, while the padded drivers' seat features width and angle adjustable arms and stylish accent stitching.



Standard Safety Features

Rear anti-tippers, front and rear reflectors and a manual brake disengagement lever are standard on Colibri scooters. The reinforced frame enhances the durability of the scooter through daily travels.



Dual Options for Tires and Batteries

Sleek double-spoked black rims on non-marking low profile tires available in 200 x 50 and 210 x 65 sizes. 12 Ah and 18 Ah battery options provide convenient charging on and off the scooter and a travel range of up to 16 km.



Five Lightweight Components

- Seat 22 lbs
- Front 3W 29 lbs
- Front 4W 34 lbs
- Rear 21 lbs
- Battery Pack 12 Ah 22 lbs
- Battery Pack 18 Ah 32 lbs
- Basket 1.45 lbs



Leo™ 3-Wheel & 4-Wheel Indoor/Outdoor

Model Numbers:

	3-Wheel	4-Wheel
Blue	LEO-3B	LEO-4B
Red	LEO-3R	LEO-4R
Silver		LEO-4S

Available in three and four-wheel platforms, the Leo scooter is the ideal indoor/outdoor scooter due to its compact size and long-range capacity. Engineered with the highest regards to safety, function and quality, the Leo scooter is the optimal choice for those who value their independence. Features such as delta-steering handles, an integrated lighting system and a comfort-adjustable tiller and seat ensure that any journey on the Leo scooter is comfortable and safe. Includes 36 Ah batteries.

Features

- Standard with 36 Ah batteries providing up to a 36 km range.
- Tool-free disassembly into five lightweight parts for travel and transportation.
- Compact size ensures easy control and maneuverability.
- Built-in splash guards protect the motor and electronics.
- Black, pneumatic tires are non-scuffing and provide a soft ride.
- 10" black pneumatic tires with polished aluminum rims.
- Two-step disengaging lever guards against unintentional disengagement out of drive mode.
- 360° rotating seat is height adjustable and features width and angle adjustable arms and stylish accent stitching.
- Full lighting system. Efficient LED lights increase visibility and improve scooter range.
- Provider-programmable electronics to reverse the throttle direction and adjust driving performance.
- Infinitely adjustable tiller with ergonomic delta-steering handles.
- Easy to use, lit control panel with programmable audible and visual signals.
- Dual rear accessory ports accommodate the use of multiple accessories at one time.
- 350 lbs. weight capacity.

Leo™ Scooter Features



Infinite Tiller Adjustment

The angle of the tiller easily adjusts to ensure a comfortable driving position and to allow for easy control panel interaction and viewing.



Comfort Seating

Height adjustable seat features width and angle adjustable arms, fore/aft adjustment and 360° of rotation. Dual accessory ports allow for the use of two rear accessories, if required.



Brake Light/Indicators

Rear brake lights engage upon the release of the throttle. Left and right direction indicators automatically stop blinking after 30 seconds.



Intuitive Control Panel

Designed to ensure easy and accurate operation with domed buttons, LED lights and audible signals as well as a lit battery indicator showing remaining battery power.



Splash Guards

Built-in splash guards protect the motor and electronics and help prevent unnecessary servicing.



PEGASUS3M-R

PEGASUS4M-W

Pegasus^{METRO} 3-Wheel & 4-Wheel Outdoor

Model Numbers:

	3-Wheel	4-Wheel
Blue	PEGASUS3M-B	PEGASUS4M-B
Red	PEGASUS3M-R	PEGASUS4M-R
Silver	PEGASUS3M-S	PEGASUS4M-S
White	PEGASUS3M-W	PEGASUS4M-W
Green	PEGASUS3M-G	PEGASUS4M-G

The Pegasus^{METRO} three and four-wheel scooters are equipped with features designed to make your outdoor driving experience an enjoyable one. Proprietary Speed Reduction Technology ensures safe maneuvering around corners, and an intuitive control panel with visual and audible signals combine with programmable electronics to make the Pegasus^{METRO} exceptionally easy to drive.

Features

- Invacare’s Speed Reduction Technology (SRT).
- Ergonomic design to encourage better alignment of the head, shoulders and elbows for a more relaxed driving position.
- Infinite adjustment of the tiller between 70° and 90°.
- 4-pole motor provides reliable performance and a 12 km/h top speed.
- Full front and rear lighting system with automatic brake lights, turning signals and super bright, dual bulb headlights.
- Intuitive control panel with domed buttons, LED indicators and programmable audible signals for concise, easy operation.
- Contoured captain’s chair is height adjustable and offers enhanced multi-position recline, fore/aft adjustment and armrest height, width and angle adjustment capabilities.
- Programmable electronics.
- Adjustable rear and front spring suspension provide smooth travel on uneven terrain.
- 11” black pneumatic tires with polished aluminum rims.
- Two-step brake release lever prevents accidental disengagement out of drive mode.
- Motor and electronics that are protected from corrosion and water, to safeguard the vital workings of the scooter.
- Dual rear accessory ports.
- 300 lbs. weight capacity.



PEGASUS4P-R



Pegasus^{PRO} 4-Wheel Outdoor

Model Numbers:

	Classic Colours		Premium Colours
Blue	PEGASUS4P-B	Electric Blue	PEGASUS4P-EB
Red	PEGASUS4P-R	Carmin Red	PEGASUS4P-CR
Silver	PEGASUS4P-S	Iconic Chrome	PEGASUS4P-IC
White	PEGASUS4P-W	Luxor Beige	PEGASUS4P-BG
Green	PEGASUS4P-G	Crystal Copper	PEGASUS4P-CC

The Pegasus^{PRO} is a powerful upgrade to the Pegasus^{METRO} and delivers a luxury trim level for discerning riders looking to elevate their driving experience. Equipped with a backlit LCD display, the Pegasus^{PRO} provides constant feedback about the performance of the scooter while driving. Additionally, an enhanced battery capacity, larger wheels and new double wishbone suspension combined, make the Pegasus^{PRO} the ideal solution for those who want superior performance and functionality from their scooter.

Features

- An interactive LCD screen – The scooters' information center displays the time of day, temperature, battery gauge, system messages, odometer and speedometer.
- Programmable electronics, adjustable audible signal volumes and customizable LCD screen brightness.
- Invacare's Speed Reduction Technology (SRT).
- Ergonomic design to encourage better alignment of the head, shoulders and elbows for a more relaxed driving position.
- Infinite adjustment of the tiller between 70° and 90°.
- 4-pole motor provides reliable performance and a 15 km/h top speed.
- Full front and rear lighting system with automatic brake lights, turning signals and super bright, dual bulb headlights.
- Contoured captain's chair is height adjustable and offers enhanced multi-position recline, fore/aft adjustment and armrest height, width and angle adjustment capabilities.
- Vertical articulating double wishbone suspension which delivers superior control, performance and comfort.
- 12" black pneumatic tires with polished aluminum rims.
- Two-step brake release lever prevents accidental disengagement out of drive mode.
- Motor and electronics that are protected from corrosion and water, to safeguard the vital workings of the scooter.
- 350 lbs. weight capacity.

STYLE IS FOREVER. INVACARE SCOOTERS.



Comet^{PRO} 4-Wheel Outdoor

Model Numbers:

	Classic Colours		Premium Colours
Blue	COMET4P-B	Electric Blue	COMET4P-EB
Red	COMET4P-R	Carmine Red	COMET4P-CR
Silver	COMET4P-S	Iconic Chrome	COMET4P-IC
White	COMET4P-W	Luxor Beige	COMET4P-BG
Green	COMET4P-G	Crystal Copper	COMET4P-CC

Comet^{PRO} four-wheel scooters offer power and performance in a premium package. Combining large 13" pneumatic tires with superior front wishbone suspension and central rear spring suspension, the Comet^{PRO} is ideal for tackling outdoor terrain. The Comet^{PRO} also features an interactive and informative LCD screen, displaying your speed, battery status, time, temperature and system messages.

Features

- An interactive LCD screen – The scooters' information center displays the time of day, temperature, battery gauge, system messages, odometer and speedometer.
- Programmable electronics, adjustable audible signal volumes and customizable LCD screen brightness.
- Invacare's Speed Reduction Technology (SRT).
- Ergonomic design to encourage better alignment of the head, shoulders and elbows for a more relaxed driving position.
- Infinite adjustment of the tiller between 70° and 90°.
- 4-pole motor provides reliable performance and a 15 km/h top speed.
- Full front and rear lighting system with automatic brake lights, turning signals and super bright, dual bulb headlights.
- Contoured captain's chair is height adjustable and offers enhanced multi-position recline, fore/aft adjustment and armrest height, width and angle adjustment capabilities.
- Vertical articulating double wishbone front suspension which delivers superior control, performance and comfort.
- 13" black pneumatic tires with polished aluminum rims.
- Two-step brake release lever prevents accidental disengagement out of drive mode.
- Motor and electronics that are protected from corrosion and water, to safeguard the vital workings of the scooter.
- 350 lbs. weight capacity.

STYLE IS FOREVER. INVACARE SCOOTERS.



COMET4HDU-G

Comet^{ULTRA} 4-Wheel Outdoor

Model Numbers:

	Classic Colours	Premium Colour Kits	
Blue	COMETHDU-B	Electric Blue	1637418
Red	COMETHDU-R	Carmine Red	1637417
Silver	COMETHDU-S	Iconic Chrome	1637419
White	COMETHDU-W	Luxor Beige	1637415
Green	COMETHDU-G	Crystal Copper	1637416

The Comet^{ULTRA} was developed with full regard for comfort, stability and safety. Its robust construction supports up to 485 lbs. and features a 26" wide seat, while its unique suspension system is adjustable to reduce vibration and increase driving satisfaction. The Comet^{ULTRA} incorporates a reinforced seat base, seat column and armrests which offer reliable support when getting on and off the scooter.

Features

- 26" wide contoured captain's seat with pneumatic lumbar support that was developed in conjunction with therapists to ensure an optimal seating position and increased comfort.
- An interactive LCD screen – The scooters' information center displays the time of day, temperature, battery gauge, system messages, odometer and speedometer.
- Programmable electronics, adjustable audible signal volumes and customizable LCD screen brightness.
- Invacare's Speed Reduction Technology (SRT) .
- Ergonomic design to encourage better alignment of the head, shoulders and elbows for a more relaxed driving position.
- Infinite adjustment of the tiller between 70° and 90°.
- Powerful 4-pole motor provides a 10 km/h top speed.
- Full front and rear lighting system with automatic brake lights, turning signals and super bright, dual bulb headlights.
- Vertical articulating double wishbone front suspension which delivers superior control, performance and comfort.
- 13" black pneumatic tires with polished aluminum rims.
- Two-step brake release lever prevents accidental disengagement out of drive mode.
- Robust motor and electronics that are protected from corrosion and water, to safeguard the vital workings of the scooter.
- 485 lbs. weight capacity.

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} and Comet^{ULTRA} Scooter Features



Intuitive Control Panel

Pegasus^{METRO}

Designed to ensure easy and accurate operation with domed buttons, LED lights and audible signals as well as a lit battery indicator displays remaining battery power.



LCD Control Panel

Pegasus^{PRO}, Comet^{PRO}, Comet^{ULTRA}

The LCD screen displays the speed, distance travelled, battery status, time, temperature and system messages.

The screen brightness and the volume of the panel's audible signals can be adjusted through the control panel.



Ergonomic Steering

The design of the tiller and delta handles encourages better alignment of the head, shoulders and elbows. This, combined with 70°-90° of infinite tiller adjustment places the body into a more relaxed driving position, decreasing strain and increasing comfort.



Advanced Suspension

Pegasus^{PRO}, Comet^{PRO}, Comet^{ULTRA}

Vertical articulating double wishbone front suspension delivers superior control, performance and comfort. Together with adjustable suspension in the rear frame, this advanced system provides a smooth drive, improved traction and superior stability over varying terrains.



Speed Reduction Technology

Maintain optimal control of your scooter with exclusive Speed Reduction Technology (SRT) which slows travel speed while cornering.

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} and Comet^{ULTRA} Scooter Features



Programmable Electronics

Invacare scooters are fully programmable by a provider technician to adjust driving performance, maximum speed, sleep timers and audible signals. The direction of the throttle can be reversed with a programmer.



Comfort Seating

Height adjustable seat features width and angle adjustable arms, fore/aft adjustment and 360° of rotation. The Comet^{ULTRA} seat is 26" wide and is enhanced with pneumatic lumbar support. Dual accessory ports allow for the use of two rear accessories.



Brake Light/Direction Indicators

Safety and visibility are enhanced with rear brake lights, that engage upon the release of the throttle, left and right direction indicators, which automatically stop blinking after 30 seconds.



Two-Step Brake Disengagement Lever

Prevents the accidental release out of drive mode.



Reliable Components

The motor and gearbox have been extensively tested for reliability and durability.



Secure Transport Brackets

Frame-welded transport tie-down brackets provide four points for securing your scooter during transportation.

Shroud Options



Classic Colour Shroud Options

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} & Comet^{ULTRA}



Jasper Red



Jade Green



Diamond Silver



Marble White



Sapphire Blue

Premium Colour Shroud Options

Pegasus^{PRO} & Comet^{PRO}



Crystal Copper



Iconic Chrome



Carmine Red



Electric Blue



Luxor Beige



Polygon Dusk



Polygon Sunrise



Hexagon Moguls

STYLE IS FOREVER. INVACARE SCOOTERS.

Scooter Accessories



Lockable Front Box

Mounts onto the front of the tiller to secure valuables. Features a 10 lb./4.5 kg capacity.

ACC260 (silver), ACC270 (blue), 1551638 (red),
1537710 (slate grey)

Leo, Pegasus, Comet & Comet HD

1637910 (charcoal grey)

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} & Comet^{ULTRA}



Lockable Rear Storage Box

Mounts into the rear accessory ports of the scooter for extra storage space.

1637913

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} & Comet^{ULTRA}



Extra Large Front Basket

Mounts onto the front of the tiller to provide three times the basket space versus the standard composite basket.

10.5" H x 17.75" W x 10.5" D

1637915

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} & Comet^{ULTRA}



Rear Basket

Large basket is ideal for storing extra items. 18" W x 13" D x 7" H.

ACC230

Suitable for all scooters.

Colibri requires 1148158, scooter accessory bracket.

Comet HD/Comet^{ULTRA} require 1541502, scooter extension tube.



Front and Rear Steel Bumpers

Mounts securely onto the scooter, providing additional protection for the frame and shroud.

1528275 - Front - 3W Pegasus, Pegasus^{METRO}

1638965 - Front - 4W Pegasus^{METRO}, Pegasus^{PRO}

1528298 - Front - Comet, Comet^{PRO}, Comet^{ULTRA}

1638966 - Rear - Pegasus^{METRO}, Pegasus^{PRO}

1638967 - Rear - Comet^{PRO}, Comet^{ULTRA}



Deluxe Mirror Set

Stylized front mirrors for mounting on the left and right side of the control panel. Sold as a pair.

1637911

Suitable for all scooters except the Colibri.



Standard Mirror

Removable mirror is adjustable and can be mounted on the left or right side of the control panel. Sold individually.

ACC250

Suitable for all scooters except the Colibri.



Beverage Holders

RAM Level Cup™ clamp-on drink holder.

TAGHCUP (left) Suitable for all scooters.

Foldable Cup Holder. Installs onto tiller pocket.

1521594 (right)

Suitable for all scooters except the Colibri and Leo.

Scooter Accessories



Mobile Phone Holder

Sturdy and secure bracket which adjusts to fit and hold most smartphones.

1637912

Suitable for all scooters except the Colibri.



USB Charger

Plugs into the scooters' charger port to enable the charging of one mobile device.

1637917

Pegasus^{METRO}, Pegasus^{PRO}, Comet^{PRO} & Comet^{ULTRA}



Safety Flag

6' bright orange flag makes the scooter easy to see from a distance. Comes complete with a secure mounting bracket.

ACC140

Suitable for all scooters.

Colibri requires 1148158, scooter accessory bracket.

Comet HD/Comet^{ULTRA} require 1541502, scooter extension tube.



Forearm Crutch and Cane Holder

Keeps a cane or forearm crutches stable and within reach during transit.

ACC200

Suitable for all scooters.

Colibri requires 1148158, scooter accessory bracket.

Comet HD/Comet^{ULTRA} require 1541502, scooter extension tube.



Walker Holder

Keeps a folding walker or rollator secure during travel.

ACC210

Suitable for all scooters.

Colibri requires 1148158, scooter accessory bracket.

Comet HD/Comet^{ULTRA} require 1541502, scooter extension tube.



Oxygen Tank Holder

Secure storage of up to a D-size oxygen tank for transit.

ACC220

Suitable for all scooters.

Colibri requires 1148158, scooter accessory bracket.

Comet HD/Comet^{ULTRA} require 1541502, scooter extension tube.



Windshield

Installs into the mirror ports by using either the mirrors or the hardware included.

1543160

Suitable for all scooters except the Colibri.

Scooter Accessories



Backrest Bag

Easily installs onto the scooter seat for additional storage.

PRC571720 (left)/ 1505731 (right)

Suitable for all scooters, except the Colibri and Leo.



Seat Suspension (LOW.MEDIUM.TALL)

Installs into the seat post for a softer ride.

1577212 (L)/1577213 (M)/1577214 (T) - Colibri

1541628 (L)/1541629 (M) - Leo

1640428 (L)/1640430 (M)/1640448 (T) - Comet^{PRO} to 201 lbs.

1640429 (L)/1640447 (M)/1640449 (T) - Comet^{PRO} to 285 lbs.

1640422 (L)/1640424 (M)/1640426 (T) - Pegasus^{Metro/Pro} to 201 lbs.

1640423 (L)/1640425 (M)/1640427 (T) - Pegasus^{Metro/Pro} to 285 lbs.



Scooter Cover

Heavy duty black nylon cover protects the scooter from the elements.

PRC571840 - Colibri & Leo

1518240 - Pegasus, Pegasus^{Metro/Pro}, Comet, Comet^{PRO}

Seat Cover

Protects seat from dirt rain and dirt.


1521684 (Not shown) - Pegasus, Pegasus^{Metro/Pro}, Comet, Comet^{PRO}



Enhanced Visibility Safety Vest

Highly reflective safety vest. One size.

SP1521592

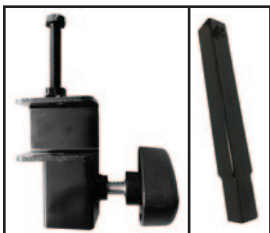

Accessory Bracket

Attaches to the seat frame of the Colibri to accommodate one rear accessory.

1148158 (left)

Extension Tube

Extends the length of the Comet HD/ Comet^{ULTRA} accessory tube to accommodate the use of rear-mounted accessories.

1541502 (right)

Product Specifications

Overall Dimensions	Colibri 3W 12Ah/ 4W 12Ah/4W 18Ah		Leo 3W/4W		Pegasus 3W/4W		Pegasus 4W		Comet 4W	
					METRO	PRO	PRO	ULTRA		
Base Width	19.25"	23.6"	26"	25.5"	26"	25.5"	26"	26"	26"	26"
Max. climbable obstacle height	1.75"/1.75"/3"	2.4"/2.5"	3"	4"	3"	4"	4"	4"	4"	4"
Length	40"	47.5"	49"/50"	52"	49"/50"	52"	57"	57"	57"	57"
Height	34.1"	36"	49.5"	49.5"	49.5"	49.5"	49.5"	50.75"	50.75"	50.75"
Usable Foot Space Length	18"/12"/12"	22"/13"	23"/13"	13"	23"/13"	13"	13.5"	13.5"	13.5"	13.5"
Seat Dimensions										
Seat-to-Floor	14.5" - 16.5"	16.5" - 19.5"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"	17.25" - 20.25"
Width	17"	19"	20"	20"	20"	20"	20"	20"	20"	26"
Depth	17.5"	16"	18.5"	18.5"	18.5"	18.5"	18.5"	18.5"	18.5"	22.5"
Back Height	15.5"	18.5"	27" - 29"	27" - 29"	27" - 29"	27" - 29"	27" - 29"	27" - 29"	27" - 29"	32" - 35"
Back Angle	N/A	N/A	90° - 130°	90° - 130°	90° - 130°	90° - 130°	90° - 130°	90° - 130°	90° - 130°	95° - 130°
Arm Dimensions										
Length	11.5"	14"	14"	14"	14"	14"	14"	14"	14"	14"
Seat-to-Top of Arms	8"	8"	7" - 10"	7" - 10"	7" - 10"	7" - 10"	7" - 10"	7" - 10"	7" - 10"	5" - 7"
Width Between Arms	16" - 19"	23" - 27"	20" - 22.5"	20" - 22.5"	20" - 22.5"	20" - 22.5"	20" - 22.5"	20" - 22.5"	20" - 22.5"	20" - 22.5"
Recommended Battery										
	2 x 12 V, 12 Ah*/	2 x 12 V, 36 Ah*	2 x 12 V, 40 Ah	2 x 12 V, 75 Ah	2 x 12 V, 40 Ah	2 x 12 V, 75 Ah	2 x 12 V, 75 Ah	2 x 12 V, 75 Ah	2 x 12 V, 75 Ah	2 x 12 V, 75 Ah
	2 x 12 V, 18 Ah*			Group 24		Group 24	Group 24	Group 24	Group 24	Group 24
Off-Board Charging System										
	2 Ah/2 Ah/5Ah	5 Ah	8 Ah	8 Ah	8 Ah	8 Ah	8 Ah	8 Ah	8 Ah	8 Ah

*included with scooter

Product Specifications

	Colibri 3W 12Ah/ 4W 12Ah/4W 18Ah	Leo 3W/4W	Pegasus 3W/4W METRO	Pegasus 4W PRO	Comet 4W PRO	Comet 4W ULTRA Pro
Speed	8 km/h (5 mph)	8 km/h (5 mph)	12 km/h (7.5 mph)	15 km/h (9.5 mph)	15 km/h (9.5 mph)	10 km/h (6 mph)
Estimated Travel Range*	11 km (7 miles) / 11 km (7 miles) / 16 km (10 miles)	36 km (22 miles)	42 km (26 miles)	55 km (34 miles)	55 km (34 miles)	58 km (36 miles)
Weight Limit (lbs.)	300**	350**	300**	350**	350**	485**
Maximum Incline Capability	6°	6°/8°	6°	10°	10°	10°
Tires						
Dimensions	200 x 50 / 200 x 50 / 210 x 65	9" x 3.5"	11" x 3.5" 4.10/3.50-5	12" x 3.5" 3.50-5	13" x 4" 4.00-5	13" x 4" 4.00-5
Tire Pressure	N/A	50 psi	50 psi	50 psi	50 psi	50 psi
Weight (lbs.)	91.6 / 97.4 / 108 (inc. batteries)	174/184 (inc. batteries)	114/125	130	138	155

Warranty Two-year limited warranty against manufacturing defects (excluding upholstered materials, padded materials and tires/wheels). Six months on batteries purchased from Invacare®. Manufacturer's warranty does not cover normal wear.

*Travel range is an estimate only based on ISO 7176-2 standards determined by testing using a 100 kg/220 lb. dummy. Travel range will vary depending on user weight, speed, battery type, battery condition and charge level, weather, terrain and condition of the drive wheels.

**Weight limit is the maximum total weight of the rider, accessories and additional items the rider may be transporting. Example: If the weight limit of the scooter is 350 lbs. and additional items equal 25 lbs., the maximum rider weight can be 325 lbs.



Yes, you can.®

Invacare Canada LP

www.invacare.ca

Canada

570 Matheson Blvd. E.

Unit 8

Mississauga, Ontario

L4Z 4G4 Canada

(800) 668-5324

All rights reserved.
Trademarks are identified
by the symbols TM, SM and [®].
All trademarks are owned
by or licensed to Invacare
Corporation unless
otherwise noted.

Specifications are subject to
change without notification.

© 2018

Invacare Canada

Form No. 08-621C

Rev 06/18