

matrx[®]
SEATING SERIES

INVACARE[®] MATRX[®]

LIBRA[®]
CUSHION



The **Invacare[®] Matrx[®] Libra[®] Cushion** — The complete solution is designed with the highest level of skin protection, positioning and adjustability for unique user needs. It is specially contoured to provide balance and stability to the pelvis, optimizing body function. Recommended for all users including those at highest risk of skin breakdown.

Outer Cover – Moisture resistant, high stretch, breathable Startex™. Hook and loop fastener strips, lifting strap, convenient zippered side pocket with easy pull zipper tab. Optional Spacetex 4000U™ cover available for extra breathability.

Inner Liner - Highly flexible Platilon liner provides an extra layer of moisture protection.

Fluid Sac – Dual layer, sealed compartments for maximum protection around boney prominences. Non-temperature sensitive, no kneading or maintenance required.

Ab/adductor contour promotes comfortable midline

leg position


Foam Base – Anatomically contoured to optimize pelvic support and positioning. High resilient foam made with Ultra-Fresh™ anti-microbial and odor protection.

Pre-ischial contour prevents sliding to reduce shear

Trochanteric shelf promotes level pelvis, lateral stability

Extended rear cushion radius helps prevent posterior pelvic tilt

Waterfall front edge accommodates tight hamstrings, foot propulsion

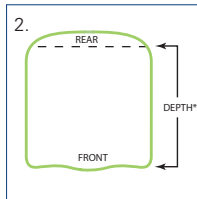

1. Modification

Leg length discrepancy cuts and depth reductions available, outer cover tailored to fit.

2. Actual cushion depth is 1 1/2" longer than stated depth due to rear cushion radius.

Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Motion Concepts

www.motionconcepts.com

USA

700 Ensminger Rd., Suite 112
Tonawanda, NY 14150
(888) 433-6818
Fax: (888) 433-6834


Canada

84 Citation Dr., Unit 1-7
Concord, ON L4K 3C1
(866) 748-7943
(905)-695-0134
Fax: 905-695-0138

© 2015

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Matrix is a registered trademark of Motion Concepts. Startex is a registered trademark of Star-Tex Corporation. Spacetex 4000U is a trademark of FR Systems International. Ultra-Fresh is a trademark of Thomson Research Associates. Specifications are subject to change without notification.
Form No. 14-399
141144

Libra Fit Kit Accessories: Designed to enhance immersion, correct or accommodate orthopedic issues.

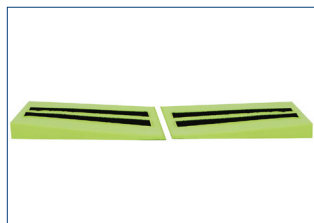


Lateral Wedges

Creates extra depth of contour for enhanced ischial relief

Tapered Wedges

Enhance adductor support, address pelvic rotation/obliquity

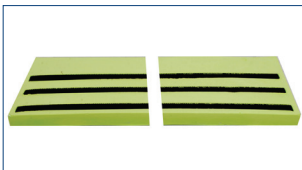

Anterior Seat Wedges

Enhances pre-ischial contour to prevent sliding

Hemi Wedge

For accommodation of hip flexion limitations, foot propulsion




Fluid Supplement

Unique strap attachment for targeted application

Abductor

To promote midline leg position

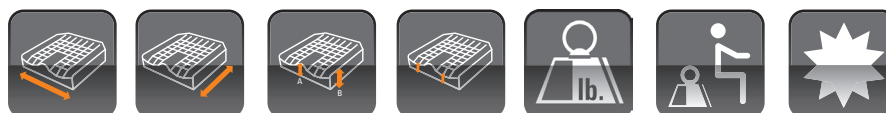

| Model # | Width | Depth | Height at Abductor/Adductor | Height at Leg Trough | Product Weight | Weight Capacity | Limited Warranty |
|-----------------|-----------|-----------|-----------------------------|----------------------|----------------|-----------------|------------------|
| LC (Standard) | 14" - 20" | 14" - 20" | 3.35" | 2.35" | 4.0 lb. and up | 300 lb. | 3 Years |
| LC (Heavy-Duty) | 21" - 24" | 18" - 22" | 3.60" | 2.60" | 4.0 lb. and up | 500 lb. | 3 Years |

